

HOW TO USE YOUR FIRST HEARING AID


- Any hearing aid cannot return level of your hearing to normal or halt further hearing deterioration, but proper use of your hearing aid can help you hear what is going on around you and let you get more enjoyment of life.
- If this is your first hearing aid, you - like most new users - will probably be surprised at the loudness of the sounds, and after using it for a short while, you may even feel that it is too noisy. The following three tips are important in overcoming your concerns.
 - 1) Read this manual and become familiar with the hearing aid's various mechanisms and how to use them correctly.
 - 2) First, use your hearing aid in a quiet place and listen to quiet sounds and give yourself sufficient time to get accustomed to the sound of your hearing aid.
 - 3) If any physical problems develop, consult your doctor.

FOR KEEPING YOUR HEARING AID LONGER

- Always keep the hearing aid clean.
Clean the hearing aid with dry cloth periodically and also remove earwax plugged in the sound outlet using eyebrush or toothbrush.
- Do not expose the hearing aid to the heat from stoves, etc. Also, avoid leaving it in places with high humidity.
- If the hearing aid is dropped in water, wipe with dry cloth and ask the dispenser for check.
- Do not pick the microphone with nail, pin, etc., because this may damage the microphone.
- Do not disassemble or attempt any repairs by yourself. Ask the dispenser for repair.

OPERATING INSTRUCTIONS


Parts and Controls

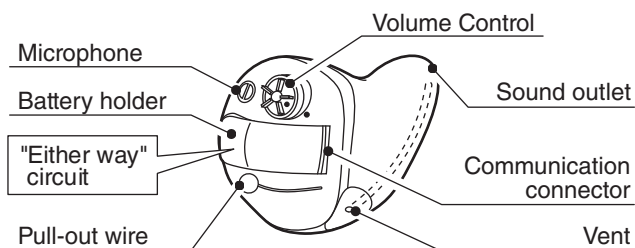

- Depending on the shape of the user's ear, there may be no vent, or the vent may be in a different position.
- * The communication connector serves for adjustment of the hearing aid by the hearing aid professional.

Battery Placement

1. Peel off the sticker of the battery (PR41 or 312 zinc-air battery) before loading.
2. Open the battery holder.
3. Insert a PR41 or 312 battery sideways into the battery holder. **You do not have to worry about the polarity.**
4. Close the battery holder.

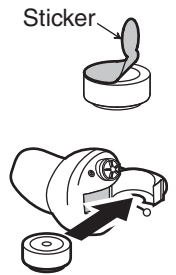


Battery life (PR41 or 312, for continuous use)

Approx. 150 hours

- Battery life will change depending on the operating condition.

When remove the battery, grasp it and pull it towards you. Be sure to remove the battery when the hearing aid is not to be used for long time.

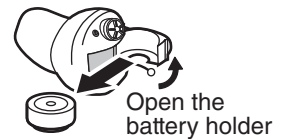


Inserting Hearing Aid in Your Ear


This hearing aid does not have a power switch. The battery holder also functions as power switch.

- Power/OFF: Open the battery holder.
Power/ON: Close the battery holder firmly.


1. Open the battery holder slightly, so that the hearing aid is set to the OFF condition.

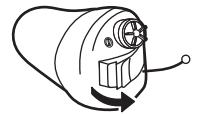

2. Hold the hearing aid so that the pull-out wire is at the rear, lightly pull the auricle back and insert the hearing aid straight into the ear.

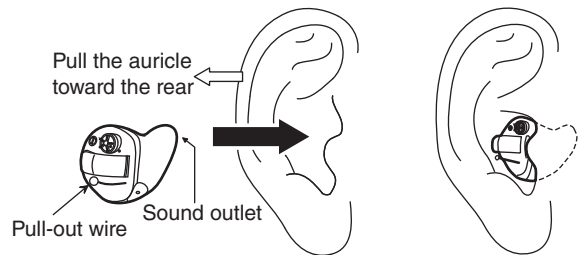


- If the hearing aid is hard to insert in a straight motion, slowly rotate it in the direction indicated by the arrow in the illustration, and lightly push it in.

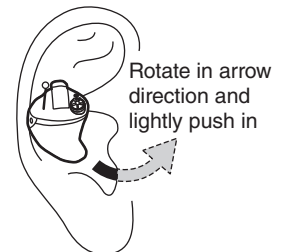

Note: If the hearing aid is inserted with the battery holder closed (ON), acoustic feedback may occur.

3. When the hearing aid is properly seated in your ear, close the battery holder to turn it on.

4. Adjust the volume control to a suitable level. Turn the control forward to increase the volume, or backward to decrease the volume.

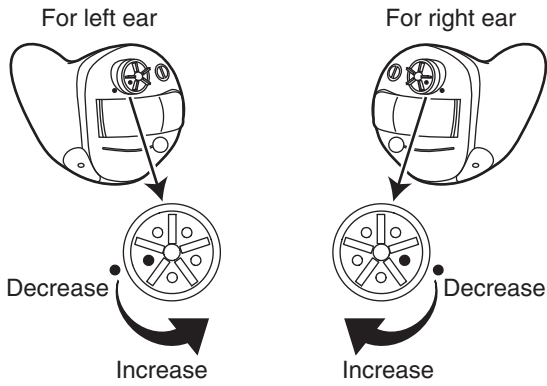

- Turn the volume control while lightly pushing it with the ball of your finger.

Removing Hearing Aid

Open the battery holder (power off), then carefully pull the wire to remove the hearing aid (if not, whistling may result). Pull the wire first sideways and then up.

TECHNICAL DATA (According to ANSI standard S3.22 2003)

Maximum OSPL90	117 dB (3000 Hz)
HFA-OSPL90	114 dB
HFA-full-on Acoustic Gain	45 dB
Reference Test Gain	37 dB
Frequency Range	200 Hz to 5000 Hz
Equivalent Input Noise Level	18 dB
Total Harmonic Distortion	500 Hz: 1% 800 Hz: 1% 1600 Hz: 1%
AGC (Attack / Recovery Time)	Input AGC Attack: 8 msec or 256 msec Recovery: 8 msec or 256 msec
	Output AGC Attack: 0.25 msec Recovery: 0.25 msec
Battery Type/Supply Voltage	312/1.4 V
Battery Current/Battery Life	0.95 mA/150 hours

(Typical value)

Turn-on precaution


If the hearing aid is turned on again immediately after turning it off, there may be no sound, but this is not a defect. In such a case, turn the power off (open battery holder) and wait 10 seconds. Then turn the power on again (close battery holder).

Adjusting the Controls

This hearing aid is a digital type. Any adjustments must be done at the dealer. For details, please contact the dealer.

Cleaning Your Hearing Aid

We recommend that you clean your hearing aid twice a week (more or less needed) to remove any earwax that may build up as a result of regular use. Too much earwax build up can prevent your hearing aid from working. Hold the hearing aid as the illustration, brush off the earwax of the sound outlet.

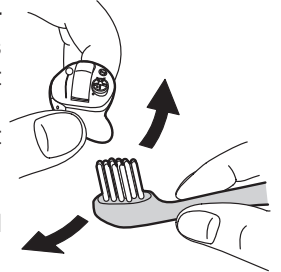


TECHNICAL DATA (According to IEC standard Pub. 60118-0-1983 Amendment 1-1994)


Reference Test Frequency	1600 Hz
OSPL ₉₀	120 dB
	500 Hz 119 dB Peak 127 dB
Full-on Acoustic Gain	50 dB
Equivalent Input Noise Level	18 dB
Total Harmonic Distortion	500 Hz: 1% 800 Hz: 1% 1600 Hz: 1%
AGC (Attack / Recovery Time)	Input AGC Attack: 8 msec or 256 msec Recovery: 8 msec or 256 msec
	Output AGC Attack: 0.25 msec Recovery: 0.25 msec
Battery Type/Supply Voltage	PR41 (312)/1.4 V
Battery Current/Battery Life	0.95 mA/150 hours

(Typical value)


OSPL90 curve

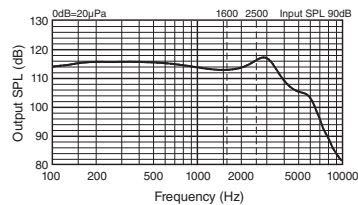


Full-on gain curve

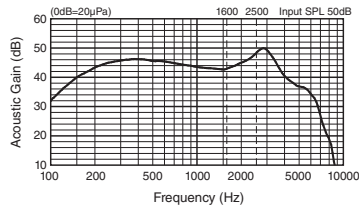


OSPL90 curve

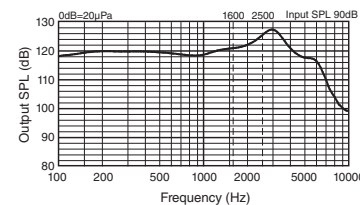

Full-on acoustic gain frequency response curve

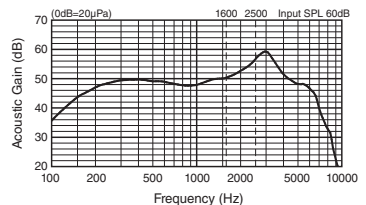

Frequency response curve

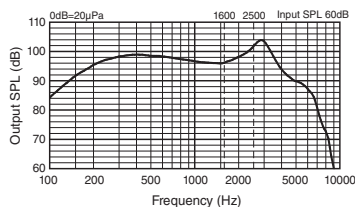

Basic frequency response curve

