

All in one model featuring superb sound quality,
Noise Reduction, Feedback Canceller and
Pulse Noise Suppressor...
All for an affordable price.

HM-06

**Digital Hearing Aid with
trimmer type volume control**

Simple

Ready-to-wear design ensures an easy fit.
Simply insert into either ear.
The device power will power off automatically
when returned to the carrying case.

Significantly improved sound in any location


Features include Noise Reduction,
Feedback Canceller and Pulse Noise
Suppressor.

Ready-to-wear

Can be used from the day of purchase.

Affordable

Superior quality at an affordable price.
Designed to be worn in either ear.

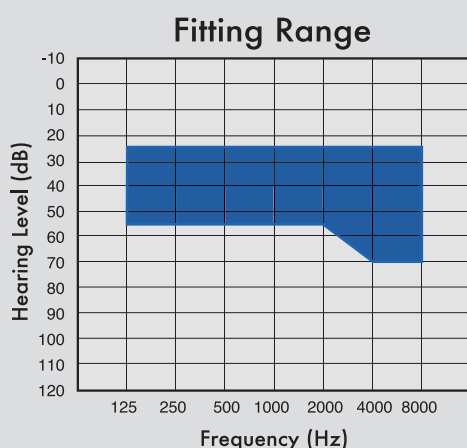


Specifications

TECHNICAL DATA (According to IEC 60118-7:2005, ANSI/ASA S3.22-2014)

OSPL90	Peak	107dB
	HFA	106dB
Full-on gain	Peak	27dB
	HFA	22dB
Reference test gain		22dB
Bandwidth (Frequency range)		$f_1 < 200\text{Hz}$ $f_2 > 5000\text{Hz}$
Equivalent input noise		18dB*
Total harmonic distortion	500Hz	1%
	800Hz	1%
	1600Hz	1%
Battery current		0.95mA
Battery type		PR41(312)
Dimensions / Weight		13.5 x 24.9 x 10.3mm / 1.6g (excluding battery)


* Low-Level expansion is active.

(Typical value)

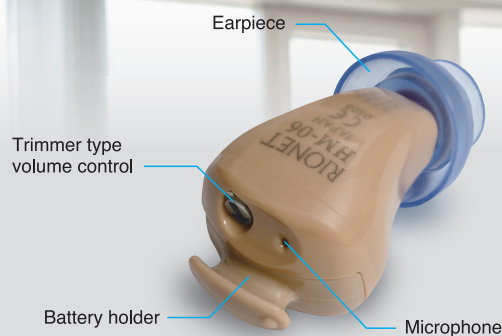


HM-06
Digital Hearing Aid with trimmer type volume control


● OSPL90 response curve

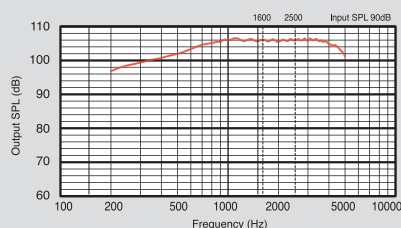

● Full-on gain response curve

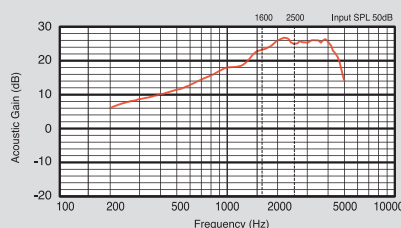

● Basic frequency response curve

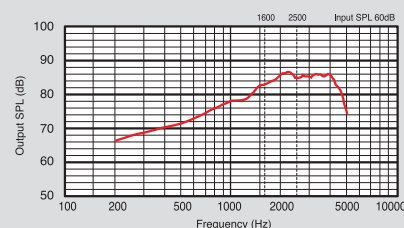

Storage of the hearing aid

To store your hearing aid, simply insert it into the case as illustrated below.
There is no need to remove the battery or open the battery holder.
Once the case is closed, the power is automatically turned off.



Specifications subject to change without notice

Made in Japan

RION CO., LTD.
URL: <http://www.rion.co.jp/english>

3-20-41 Higashimotomachi, Kokubunji, Tokyo 185-8533, Japan
Telephone: +81-42-359-7862 Fax: +81-42-359-7441

Distributed by: